ALyS series

ALyS™ 505N

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
with IL-2 1000 IU/mL

ALyS505NK
-AC1000

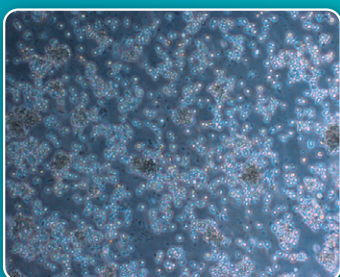- Specially designed for human NK cells
- Achieves much number of NK cells and a high NK cell ratio

Serum-free medium for NK cells

ALyS™ 505NK-AC&EX

ALyS505NK-AC and -EX are serum-free medium optimal for NK cell therapy.

- Achieve a highly-pure NK cell rate.
- It contains human derived proteins.
Any other animal derived component free.

Performance data of ALyS505NK-AC&EX *

* ALyS505NK-AC1000 / ALyS505NK-EX1000

Fig.1 Robust proliferation

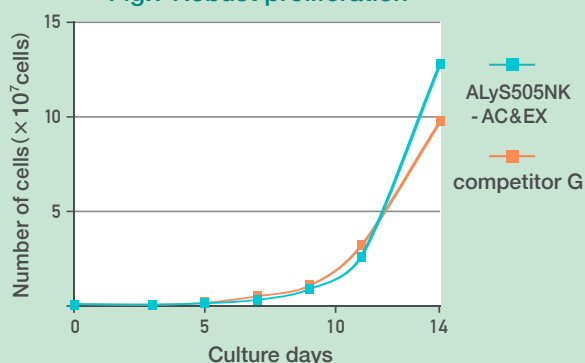
Fig.2 Number of NK cells at the day 14th

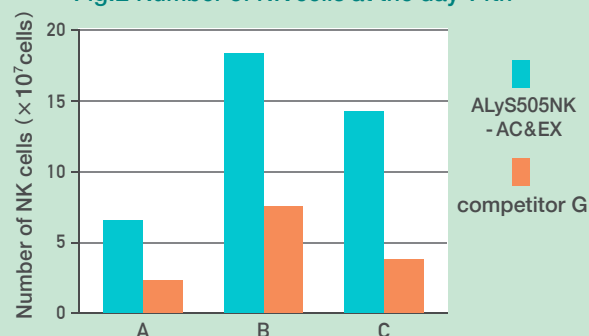
Fig.3 Ratio of NK cells and T cells


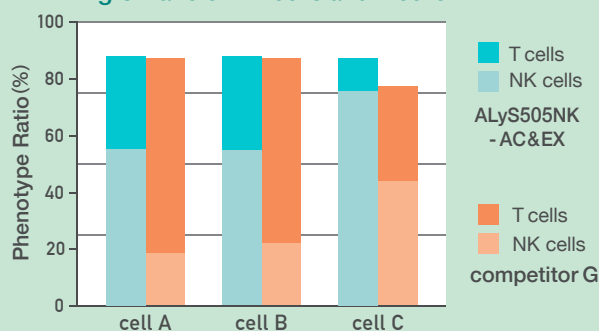
Fig.4 Cytotoxic effect of NK cells

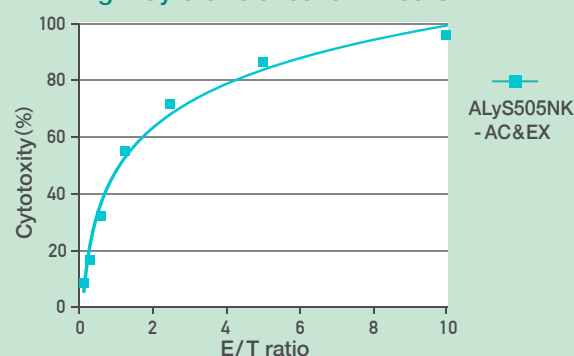Code (CSTI / NIPRO)	Product Name	Description	Volume
01600P02 / 87-554	ALyS505NK-AC	NK cell culture medium [for activation] (without IL2)	200mL (P)
01610P02 / 87-555	ALyS505NK-AC1000	NK cell culture medium [for activation] (with IL2 1000 IU/mL)	200mL (P)
01400P10 / 87-556	ALyS505NK-EX	NK cell culture medium [for expansion] (without IL2)	1,000mL (P)
01400C10 / 87-558			1,000mL (CB)
01410P10 / 87-557	ALyS505NK-EX1000	NK cell culture medium [for expansion] (with IL2 1000IU/mL)	1,000mL (P)
01410C10 / 87-559			1,000mL (CB)

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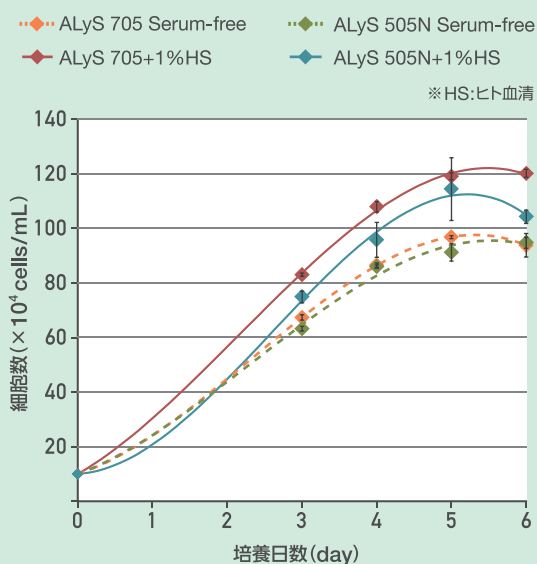
CiMS™

for Human Mesenchymal Stem Cells (hMSCs)

Basal Medium


CiMS-BM

Xeno-free
Supplement

CiMS-sXF

Animal-free
Supplement

CiMS-sAF

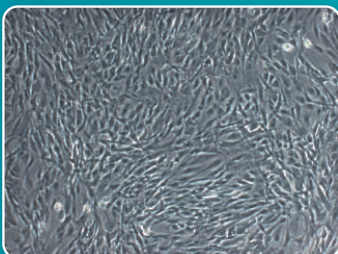- Specially designed for hMSCs
- No pre-coating required
- Cell culture tested by hMSCs

CiMS™-BM


CiMS™-sXF CiMS™-sAF

- CiMS is specially formulated culture medium to proliferate various type of hMSCs in an undifferentiated state without pre-coating of culture vessels.
- CiMS is comprised of the basal medium (CiMS-BM) and the supplement (CiMS-sXF or CiMS-sAF). Please select either Xeno-free (CiMS-sXF) or Animal-free (CiMS-sAF) supplement.
- hMSCs grown in serum-supplemented media can be transitioned into CiMS directly without adaptation.


Xeno-free: It contains human derived component. Any other animal derived component-free.
Animal-free: animal derived component-free.


1. Cell morphology (hMSCs from adipose tissue)

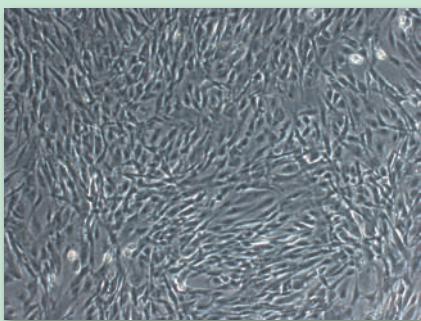
2. Robust proliferation

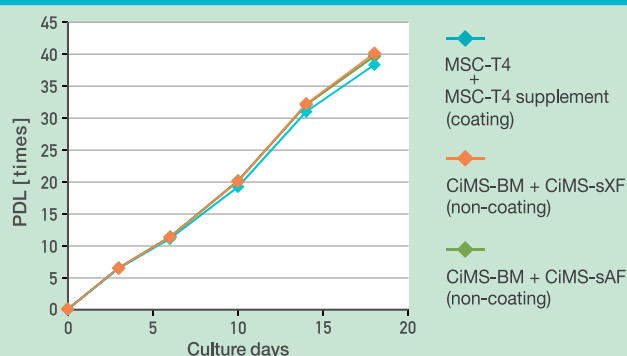
3. Multiple lineage differentiation potential

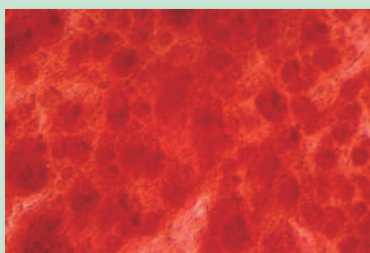Alizarin Red S

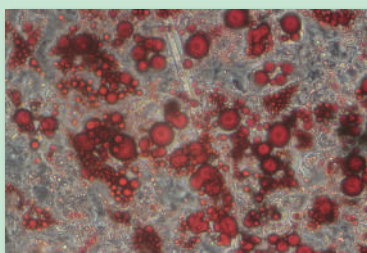Oil Red S

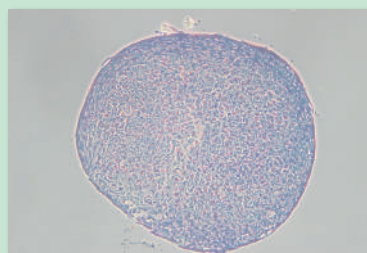Toluidine Blue

Code (CSTI / NIPRO)	Product Name	Description	Storage	Volume
A2G00P05C / 87-070	CiMS-BM	Basal Medium for hMSCs	2 – 8 °C	500mL (P)
A2G10P1CC / 87-071	CiMS-sXF	Xeno-free Supplement for CiMS-BM	below -20 °C	10 mL (P)
A2G20P1CC / 87-072	CiMS-sAF	Animal-free Supplement for CiMS-BM	below -20 °C	10 mL (P)

(P) PET bottle

※ For Research use only.

HFDM™-1


Serum-free medium for Human Fibroblast

HFDM-1(+)

rh-EGF (+)

10 ng/ml rh-EGF

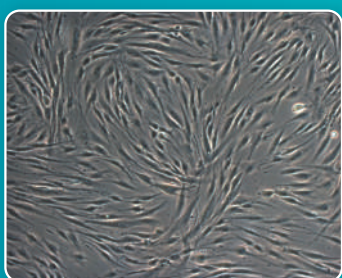- ❑ Specially designed for human normal dermal fibroblast cells
- ❑ No cell attachment factors are required
- ❑ Animal derived materials-free

HFDM™-1

- **HFDM-1** medium is a serum-free culture medium suitable for normal human fibroblast (HF) (Fig.1).
- **HFDM-1** medium contains human recombinant insulin. Human recombinant epidermal growth factor (EGF) is also included in HFDM-1(+) for the purpose of promoting cell proliferation. Ingredients derived from animal or human are not included.
- Addition of small amount of human serum (0.5-2.0%) to **HFDM-1** medium drastically enhance the growth rate of HDF cells (Fig.2).
- **HFDM-1** medium is passed with strict quality standard.

Fig.1 Human fibroblast cell growth with HFDM-1


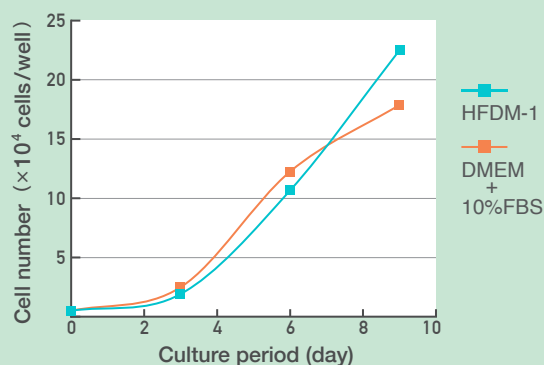
Fig.2 Additive effect of human serum

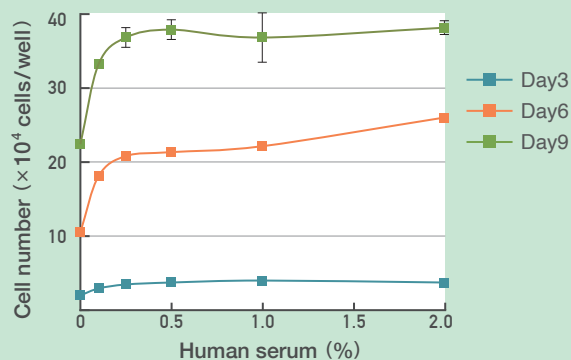Code (CSTI / NIPRO)	Product Name	Description	Volume
2102P05 / 87-160	HFDM-1 (+)	Serum-free medium for human fibroblast [with 10 ng/ml rh-EGF]	500mL (P)
2102P10 / 87-162			1,000mL (P)

(P) PET bottle

For Research use only

Artificial serum

Serum substitute

(Supplement for serum-free culture media)

Culture media
supplement

Artificial serum

Xeno-free

Culture media
supplement

Artificial serum

Animal-free




- For human Mesenchymal Stem Cells, human Fibroblasts, human Lymphocytes and other cell lines.
- A small amount of Artificial serum, 5% or less, supports cell growth as same as FBS.
- No lot-to-lot variation

Artificial serum (Xeno-free) Artificial serum (Animal-free)



- Artificial serum is developed as a substitute of serum (e.g. FBS).
- Additional "Artificial serum" to serum-free media enhance cell growth.
- Recommended concentration is 5% or less. Optimal concentration should be determined for specific use and cell type.
- Both Xeno-free and Animal-free are available.

1. Adhesion cells (Fibroblast)

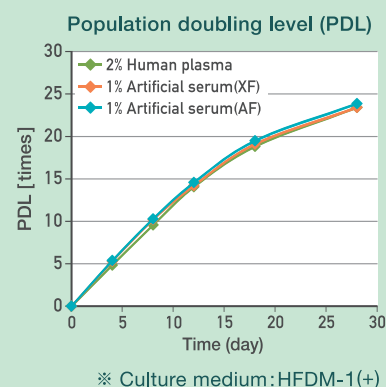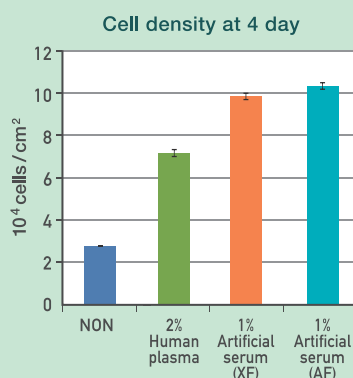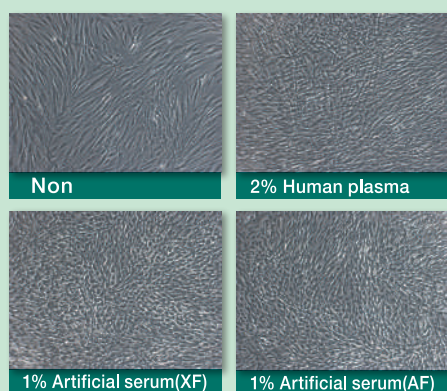2. Suspension cells (DAUDI)

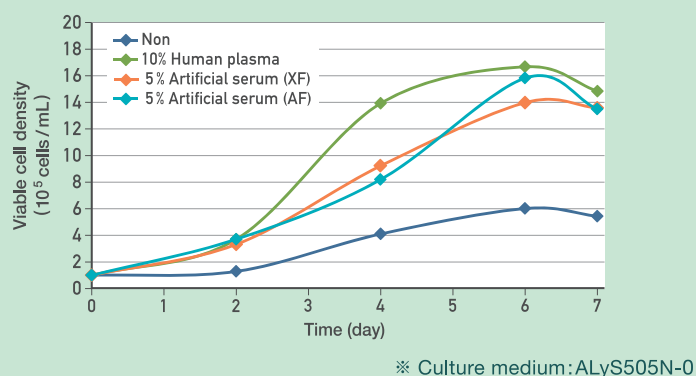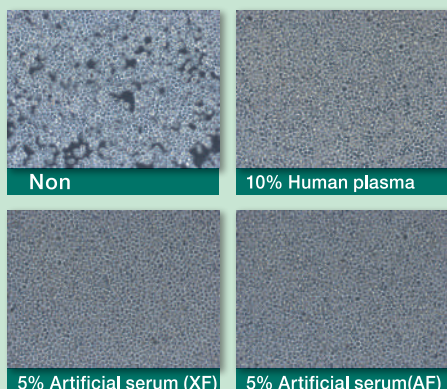Code (CSTI / NIPRO)	Product name	Description	Storage	Volume
A2G10P2CC / 87-081	Artificial serum (Xeno - free)	Xeno-free Supplement	below-20°C	20mL (P)
A2G20P2CC / 87-082	Artificial serum (Animal-free)	Animal-free Supplement	below-20°C	20mL (P)


(P) PET bottle

※ For Research use only.

HemEx[®]-Type9AØ

Mouse HSC expansion basal medium

*Requires the addition of cytokines

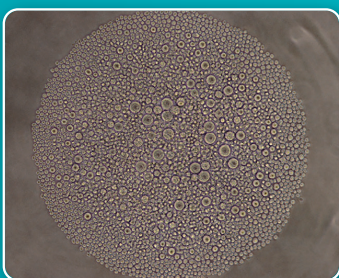- ☐ This medium is new improved version for HemEx-Type9A.
- ☐ Confirmed more mouse HSPC proliferation than conventional products.
- ☐ Supervised by Professor Hiromitsu Nakauchi and Professor Satoshi Yamazaki.

Mouse HSC expansion basal medium

HemEx®-Type9AØ *Requires the addition of cytokines

- It is possible to proliferate more undifferentiated HSC from single mouse HSC than before.
- This product developed by Dr. Yamazaki is manufactured with strict quality control in CSTI.

Mouse HSCs expansion after 10 days in culture

Fig.1 Number of cells


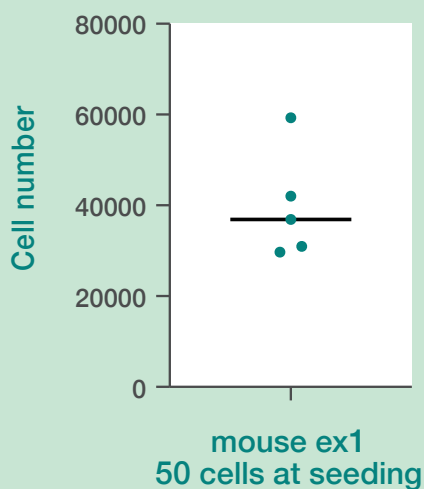
Fig.2 An example flow cytometry plot displaying mouse HSC

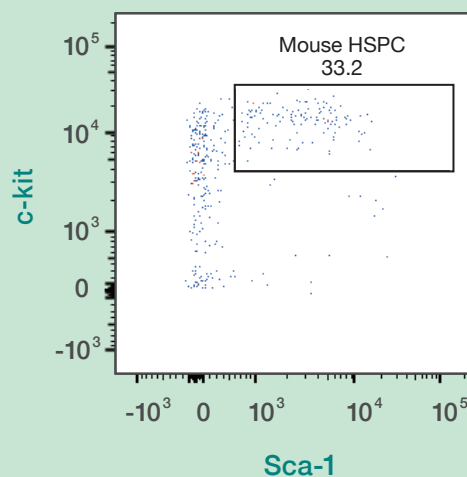Code (CSTI / NIPRO)	Product Name	Description	Volume
A5P10P01C / 87-590	HemEx-Type9AØ	Mouse HSC basal medium	100mL (P)

(P) PET bottle

※ For Research use only.

1-16-16 Seikaen, Aoba-ku, Sendai,
982-0262, Japan
URL: <http://www.cstimedia.com>
E-mail: csti-info@ml.nipro.co.jp


NIPRO CORPORATION
3-26, Senriokashinmachi, Settsu,
Osaka, 566-8510, Japan

HemEx[®]-TypeH

Human HSC expansion basal medium

*Requires the addition of chemically synthesized agonists, etc.

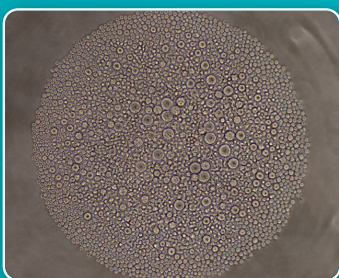- ❑ Specially designed for human HSC.
- ❑ Confirmed the proliferation of human HSPC.
- ❑ Supervised by Professor Satoshi Yamazaki.

HemEx[®]-TypeH

*Requires the addition of chemically synthesized agonists, etc.

- HemEx-TypeH is optimized for expansion of human HSC.
- This product developed by Dr. Yamazaki is manufactured with strict quality control in CSTI.
- For information on obtaining Butyzamide, a chemically synthesized agonist to be added, please contact Celaid Therapeutics Inc. (contact@celaidtx.com)
- It is necessary to add L-glutamine at a concentration of 4mM.

Human HSCs expansion after 10 days in culture

Fig.1 Number of cells


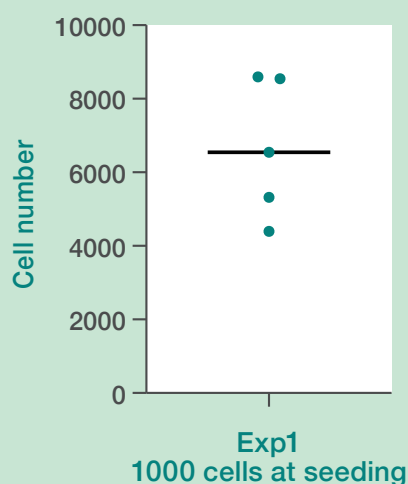
Fig.2 An example flow cytometry plot displaying human HSC

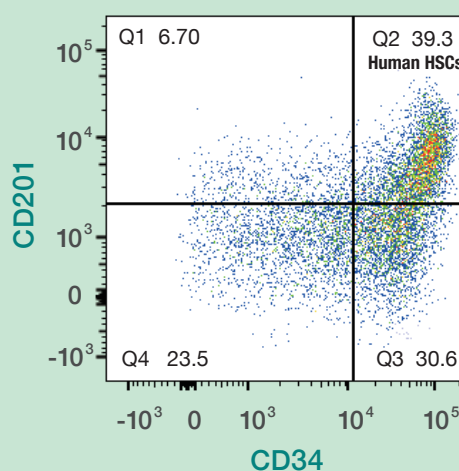Code (CSTI / NIPRO)	Product Name	Description	Volume
A2D10P01C / 87-079	HemEx-TypeH	Human HSC basal medium	100mL (P)
Supplement			
Code (CSTI / NIPRO)	Product Name	Description	Volume
1104 / 87-885	200mM L-Glutamine solution	200mM L-Glutamine solution	100mL (P)

(P) PET bottle

※ For Research use only.